

Community Face Mask

Instructions

Thank you for being part of the
mask making movement
to help our mid-level health care agencies
combat the COVID-19 pandemic.

You are making a difference, and we appreciate
you for lending your talents and time for the
cause.

For the Safety of all . . .

Please follow the guidelines listed below:

- Avoid making masks if you are sick or think you might be getting sick.
- Avoid making masks if you have been told you have a multi-drug resistant illness like MRSA or C. Diff., or if you have any open or draining wounds.
- Keep all face mask materials away from household pets, small children, or any tobacco products.
- Be sure the area in which you are working is clean (sanitize all surfaces with a disinfectant) before making masks.
- And be sure to perform hand hygiene before touching mask materials.

Thank you being a part of our goal to help our
mid-level health agencies.



SCAN THIS QR CODE WITH YOUR PHONE
TO WATCH AN INSTRUCTIONAL VIDEO

Sewing Mask provided by Canvas Innovations LLC

This is a labor of our love to our fellow Americans. These masks are **NOT** to be sold. We have put together these kits so you can help your local health care agencies. Each pack should contain 10, 50, or 100 masks. (Depending on your order)



Start by washing your hands!



1. Each mask is made of two pieces of material, plus 2 elastic straps.



2. Place the two pieces of the mask together with the darker sides facing out. (the inside color is the lighter side of fabric)



3. Sew pieces together down the center of the round side leaving about 3/8" seam allowance so you have plenty of fabric to top stitch.



4. Turn face mask pieces inside out.



5. Fold over the seam and stitch it down so the seam does not bother the nose.



6. This is what it will look like when you are done with the top stitch.




7. Next identify the larger elastic piece and sew to **TOP** edge as shown in picture.






8. Continue sewing downward and add shorter strap on bottom edge as shown in picture.

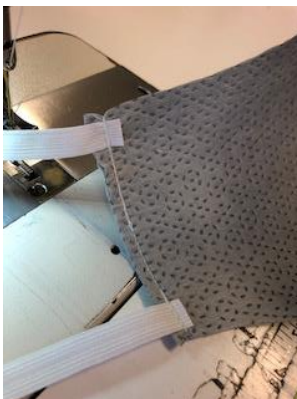

9. This is what it should look like.

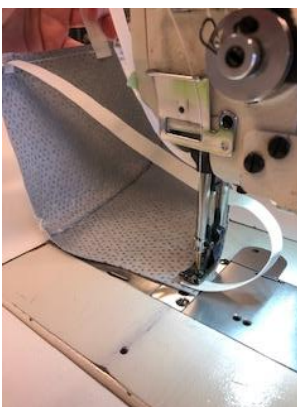

10. Now flip over and repeat on other side. Make sure to keep strap free from twists.



11. The finished mask should like this.

